GC Fuji II LC CAPSULE

RADIOPAQUE LIGHT CURED REINFORCED GLASS IONOMER RESTORATIVE IN CAPSULES

Prior to use, carefully read the instructions for use.

For use only by a dental professional in the recommended indications.

RECOMMENDED INDICATIONS

1. Restoration of Class III, V and limited Class I cavities.
2. Restoration of primary teeth.
3. Core build up.
4. Cases where a radiopaque restoration is required.
5. Geriatric applications.
6. As a base or liner.

CONTRAINDICATIONS

1. Pulp capping.
2. In rare cases the product may cause sensitivity to some persons. If such reactions are experienced, discontinue the use of the product and refer to a physician.

DIRECTIONS FOR USE

Powder / Liquid Ratio (g/g)	0,33 / 0,10
Mixing Time (sec.)	10"
Working Time (min., sec.)	3'15"
Light Curing Time (sec.)	20"
Depth of Cure (A2) (mm)	1,8

Test conditions : Temperature (23 +/-1°C)
Relative humidity (50 +/-10%)
ISO 9917-2 : 1998 (E) (Light-activated cements) (Type I)

1. CAVITY PREPARATION

- a) Prepare tooth using standard techniques. Extensive mechanical retention is not necessary. For pulp capping, use calcium hydroxide.
- b) Wash and dry but DO NOT DESICCATE. GC CAVITY CONDITIONER or GC DENTIN CONDITIONER is recommended to remove the smear layer (Fig. 1).

2. CAPSULE ACTIVATION AND MIXING

- a) Before activation, shake the capsule or tap its side on a hard surface to loosen the powder (Fig. 2).
- b) To activate the capsule, push the plunger until it is flush with main body (Fig. 3).
- c) Immediately place the capsule into a metal GC Capsule Applier and click the lever once (Fig. 4). The capsule is now activated.
Note :
The capsule should be activated just before mixing and used immediately.
- d) Immediately remove the capsule and set it into a mixer (an amalgamator) and mix for 10 seconds at high speed (+/-4,000 RPM) (Fig. 5).

3. RESTORATIVE TECHNIQUE

- a) Immediately remove the mixed capsule from the mixer and load it into the GC Capsule Applier.
- b) Make two clicks to prime the capsule then syringe (Fig. 6). The working time is 3 minute 15 seconds from start of mixing at 23°C (73.4°F). Higher temperature will shorten working time.
Note :
1) To adjust the direction of the nozzle, hold the applicator with the capsule towards you and turn the capsule body.
2) To remove the used capsule, push the applicator release button. Twist the capsule and pull upwards.
- c) Remove surface moisture but DO NOT DESICCATE.
- d) Extrude cement directly into preparation (Fig. 7). Avoid air bubbles.
- e) Form the contour and place a matrix if required.
- f) Light-cure for 20 seconds using a suitable visible light curing device (470nm wavelength) (Fig. 8). Place light source as closely as possible to the cement surface.
Note :
For cavities deeper than 1.8 mm, use a layering technique.

4. FINISHING

- Finish under water spray using superfine diamond bur, silicone point and polishing strips (Fig. 9).
Note :
If the surface has been exposed to air during curing, seal with GC Fuji VARNISH or GC Fuji COAT LC (Fig. 10).

STORAGE

Store in a cool and dark place (4-25°C / 39.2-77.0°F).
(Shelf life : 2 years from date of manufacture).

SHADE











A1, A2, A3, A3.5, A4, B2, B3, B4, C2, C4, D2 (11 shades based on Vita® shades).
Vita® is a registered trademark of Vita Zahnfabrik, Bad Säckingen, Germany.

PACKAGES

1. Single shade package : 50 capsules.
 2. Assortment package : 50 capsules (10 pcs each of A4, B2, B3, C2; 5 pcs each of B4, D2).
- * Average contents per capsule : 0.33g powder and 0.10g (0.085mL) liquid.
 - * Minimum net volume of mixed cement per capsule : 0.10mL.
- Option : GC Capsule Applier (1 piece).

CAUTION

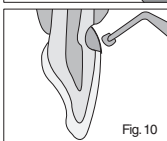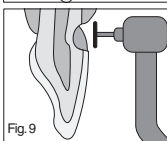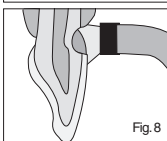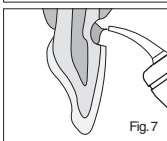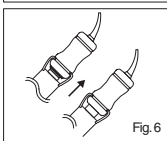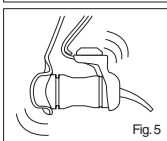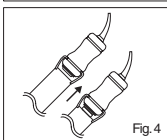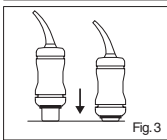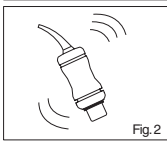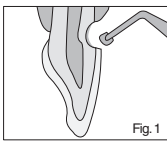1. In case of contact with oral tissue or skin, remove immediately with a sponge or cotton soaked in alcohol. Flush with water.
2. In case of contact with eyes, flush immediately with water and seek medical attention

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3737 West 127th Street, Alsip, IL 60803 U.S.A.
GC ASIA DENTAL PTE. LTD.