




COECIDE™ XL

STERILIZING AND DISINFECTING SOLUTION

PRODUCT DESCRIPTION

COECIDE XL, a 2.5% glutaraldehyde, has broad spectrum activity that is sporicidal, tuberculocidal, virucidal, bactericidal, pseudomonacidal and fungicidal. Shipped in unactivated gallons with a buffer solution for activation, **COECIDE XL** can be activated one quart at a time for economy. Activated solution can be reused for a 28-day period. **COECIDE XL** will disinfect immersed instruments at high levels in 90 minutes at 25°C. Sterilize them in 10 hours at 25°C.

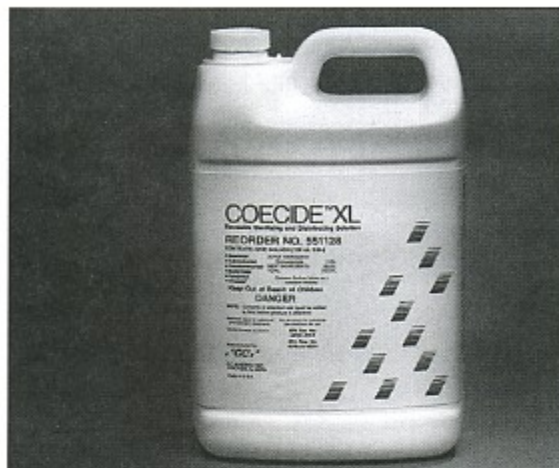

FEATURES AND BENEFITS

Stays effective for repeated use for up to 28 days after activation.

- Total gallon need not be activated—individual quarts can be poured off and activated separately.

Instruments completely sterilized in 10 hours; all bacteria, fungi, Myco-bacterium tuberculosis, even lipophilic and hydrophilic viruses are completely destroyed, including HIV-1 (AIDS Virus).

- Eliminates possibility of cross-contamination.

Contains corrosion inhibitors.

- Protects valuable instruments from harm.

Special formula has wetting and penetration ability.

- Solution reaches even the most inaccessible areas of instruments.

Packaged in marked one gallon container bottle with activator.

- Convenient and economical. . . saves you time and money.

PRODUCT CODE

551128 COECIDE XL Gallon
(3.8 L) Bottle



GC America Inc.

Advancing the Art and Science of Dentistry

© 1998 GC America Inc., 3737 West 127th St., Alsip, Illinois 60803

855000:II Rev. 9/98

<http://www.gcamerica.com>