

COMPLETE SIZE SELECTION CHART - Bold Sizes are in the Starter Kit (35) - Product Codes designate Refill Packs

Product # Producto Produto	Crown # Coronas Coroas Couronnes	Width (mm) Anchura Largura Largeur
Upper Right Central Incisor		
412064	1+.1	7.2
412065	1+.2	7.5
412066	1+.3	8.0
412067	1+.4	8.5
412068	1+.5	9.0
412069	1+.6	9.5
412070	1+.7	10.0
Upper Left Central Incisor		
412071	+1.1	7.2
412072	+1.2	7.5
412073	+1.3	8.0

Product # Producto Produto	Crown # Coronas Coroas Couronnes	Width (mm) Anchura Largura Largeur
412074	+1.4	8.5
412075	+1.5	9.0
412076	+1.6	9.5
412077	+1.7	10.0
Upper Right Lateral		
412078	2+.1	5.7
412079	2+.2	5.9
412080	2+.3	6.3
412081	2+.4	7.0
412082	2+.5	7.5
412083	2+.6	8.0
412084	2+.7	8.3

Product # Producto Produto	Crown # Coronas Coroas Couronnes	Width (mm) Anchura Largura Largeur
Upper Left Lateral		
412085	+2.1	5.7
412086	+2.2	5.9
412087	+2.3	6.3
412088	+2.4	7.0
412089	+2.5	7.5
412090	+2.6	8.0
412091	+2.7	8.3
Lower Central and Lateral		
412092	-1-2-.1	4.3
412093	-1-2-.2	4.6
412094	-1-2-.3	4.7
412095	-1-2-.4	4.9
412096	-1-2-.5	5.1
412097	-1-2-.6	5.5
412098	-1-2-.7	5.9
412099	-1-2-.8	6.2
412100	-1-2-.9	6.5
412101	-1-2-.10	6.7

Product # Producto Produto	Crown # Coronas Coroas Couronnes	Width (mm) Anchura Largura Largeur
Upper Right Cuspid		
412102	3+.1	7.2
412103	3+.2	7.4
412104	3+.3	7.6
412105	3+.4	7.7
412106	3+.5	8.0
412107	3+.6	8.4
412108	3+.7	8.9
Upper Left Cuspid		
412109	+3.1	7.2
412110	+3.2	7.4
412111	+3.3	7.6
412112	+3.4	7.7
412113	+3.5	8.0
412114	+3.6	8.4
412115	+3.7	8.9
Lower Right Cuspid		
412116	3-.1	4.9
412117	3-.2	5.4

Product # Producto Produto	Crown # Coronas Coroas Couronnes	Width (mm) Anchura Largura Largeur
412118	3-.3	5.8
412119	3-.4	6.7
412120	3-.5	7.6
Lower Left Cuspid		
412121	-3.1	4.9
412122	-3.2	5.4
412123	-3.3	5.8
412124	-3.4	6.7
412125	-3.5	7.6
Bicuspids		
412126	1	5.9
412127	2	6.0
412128	3	6.2
412129	4	6.5
412130	5	6.8
412131	6	7.2
412132	7	7.6
412133	8	8.2

2+.2	2+.3	1+.3	1+.	1+.5	1+.	1+.7
2+.	2+.5	+1.3	+1.	+1.5	+1.	+1.7
+2.2	+2.3	-1-2-.1	+3.5	+3.	3+.5	3+.
+2.	+2.5	-1-2-.7	-3.3	-3.	3-3	3-.
2	3		5		7	


STARTER PACKAGE

GC CROWNTEK™
POLYMETHYLMETHACRYLATE CROWNS


GC®
GC AMERICA INC.

GC CROWNTEK™ Technique Card


1. Select the correct size crown.

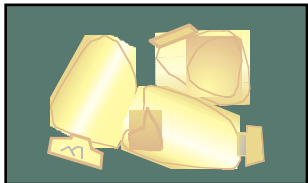


2. If needed, trim the crown margin with crown scissors or

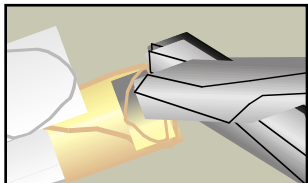

by grinding with a bur or stone (finishing the edges appropriately). Please ensure accurate seating of crown on the prepared margins.

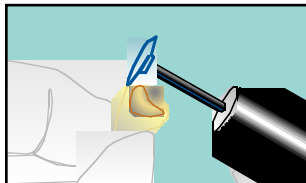


3. For good marginal adaptation, line the crown with acrylic or

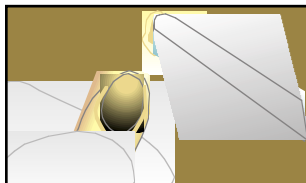


composite resin. Follow the instructions for use of resin.



If using acrylic, the crown should be removed from the preparation and resealed

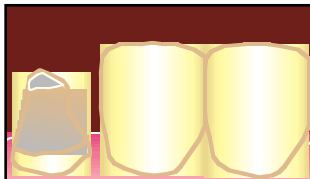

several times as the acrylic resin starts to set. (This will

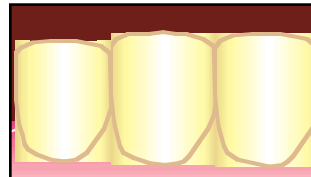


prevent locking into undercuts and help dissipate any heat build-up present during polymerization.)

4. After the resin liner has set, remove the crown from the preparation. Trim and finish margins as needed to



ensure an excellent fit. While cutting and polishing material, take care not to heat the material too much.

5. After the fit has been established, cement the provisional crown with a temporary cement. Follow the instructions for use of temporary cement.

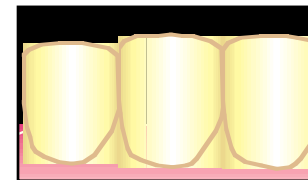

CROWNTEK Selection Chart on flipside of this card